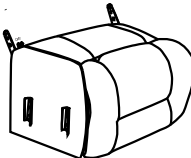
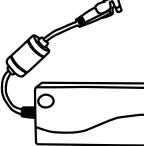
READ AND SAVE THIS INSTRUCTION FOR FUTURE USE



ASSEMBLY INSTRUCTION

PACKING DETAILS : TOTAL 2 BOXES

BOX 1/2

PART NO.	PARTS	QTY
(2)		X1
(4)		X1
(5)		X1

BOX 2/2

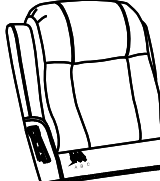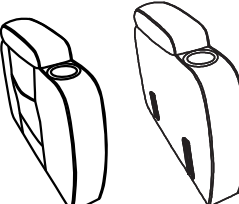PART NO.	PARTS	QTY
(1)		X1
(3) X R X L		X1 X1

We're excited for you to regain your independence with Homcom! To ensure the best experience, please read this guide carefully before using your new lift chair. Be sure to keep this manual close by, preferably in your chair pocket so that you can always refer to it as needed.

Safety Rules	4
Lift Chair Parts	6
Assembly Instructions	7
Lift and Recliner Function	8
Troubleshooting	9

SAFETY RULES

- Never set the lift chair in wet/damp areas or where dust can accumulate as these could cause an electrical fault or a mechanical dysfunction.
- Never use the lift chair inside your home if the room temperature is over 40°C.
- Never place the lift chair near a heat source (oven, chimney, furnace etc...) or in direct sunlight.

Do Not use the lift chair if:

- The lift chair is wet.
- Your hands are wet.
- Only clean the chair with a dry cloth.

Situations where you should not use the lift chair:

Switch the lift chair's system off if it gets wet or seems not to be working properly, if any part of the lift chair is damaged or malfunctioning, if any part of the lift chair appears to be damaged or wiring is exposed.

Usage:

- Indoors, avoid high temperatures and wet or damp areas.

Precaution:

- Never sit on the remote or place heavy objects on top of it.
- Do not drop the remote on the floor.
- Do not tug on the electric cable.
- Do not use the remote with wet fingers.
- Do not press or touch the buttons with sharp objects.
- Never unplug the remote without first turning the power off.
- The lift chair should rest for 30 minutes after it has been used for 30 minutes.

When the system of the lift chair is not in use:

- Make sure it is switched off.
- Unplug the lift chair's system from the wall outlet when you are not going to use it for long periods of time.
- Keep the lift chair in a clean and dry place.

Disassembling the massage armchair:


- Do not remove the back of the armchair .
- Do not try to disassemble or repair the armchair yourself.

Correct way to handle the mains plug:



- Check that you have the right voltage.
- Never use force on the cable or plug.
- Never pull the cable with wet hands.
- When removing the plug, hold the plug itself, do not pull on the cable.
- Do not use the armchair massage system if either the cable or the plug are damaged.

The Parts forming the massage armchair:


- 1) Back
- 2) Armrests
- 3) Seat
- 4) Footrest
- 5) Control unit

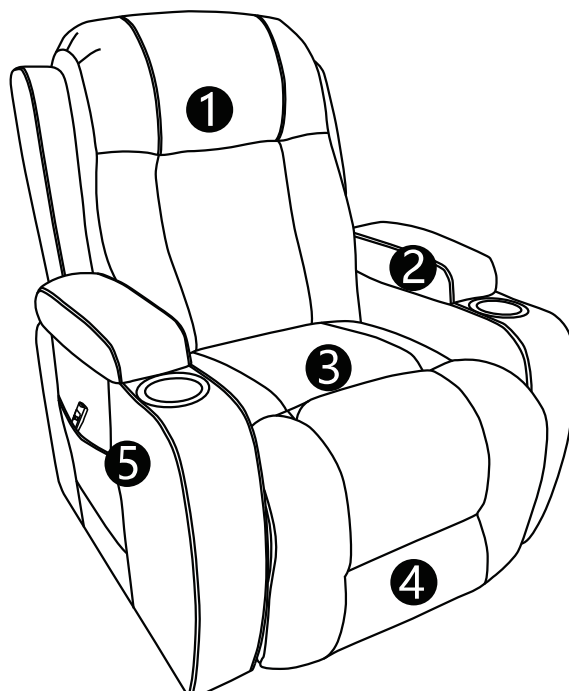
Assembly Instructions

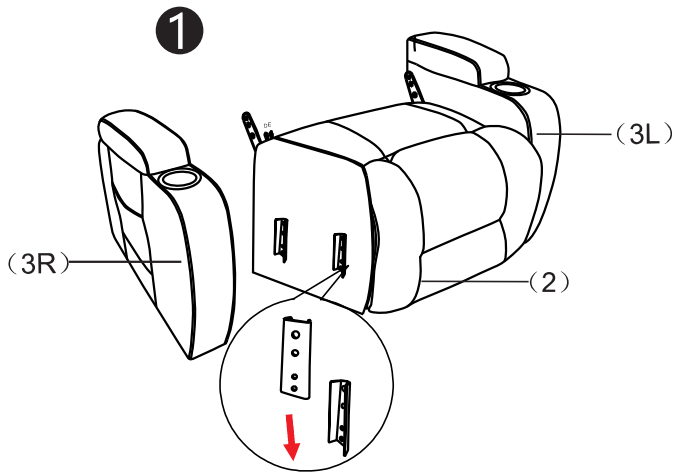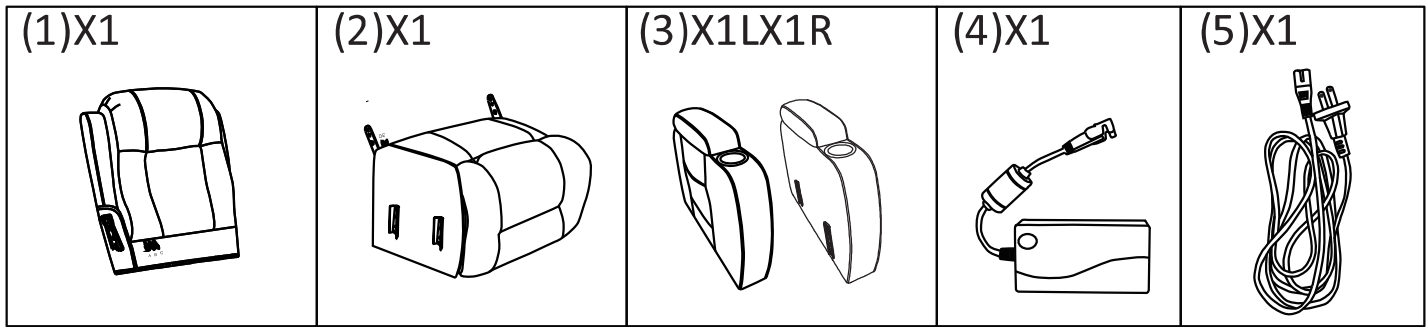
1) Insert Arms so that they fit into the seat, you will hear a “click” when they are engaged.

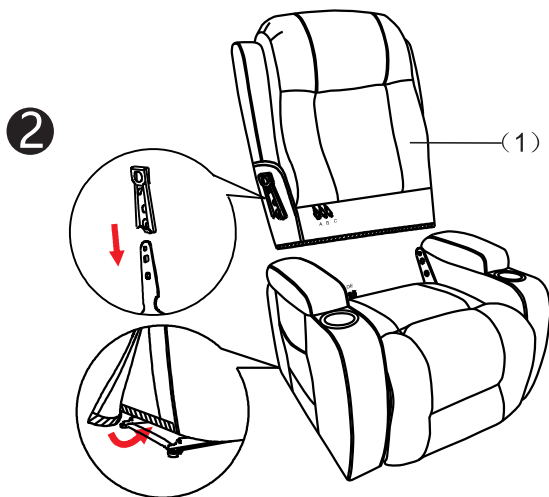2) Position Back(refer to point 1) Connect Recliner supply cable pins. Fit the pins properly in the manner shown in the image. Do not force them, they only have one correct position.

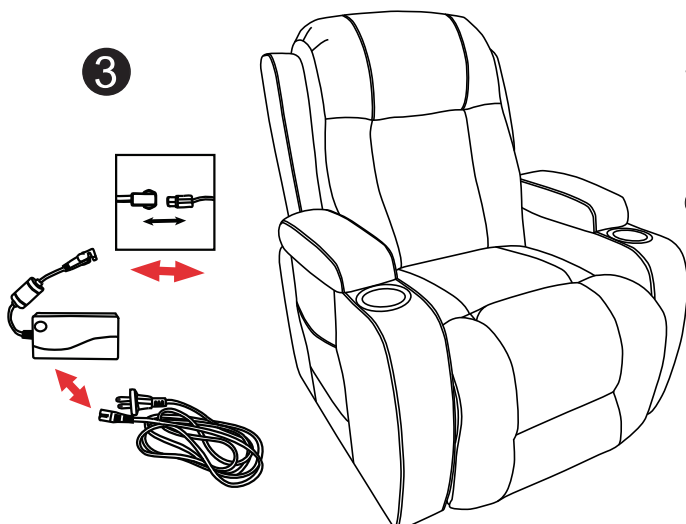3) Connect the massage and heating supply cable pins. Fit the pins properly in the manner shown in the image. Don't force them, they only have one correct position.

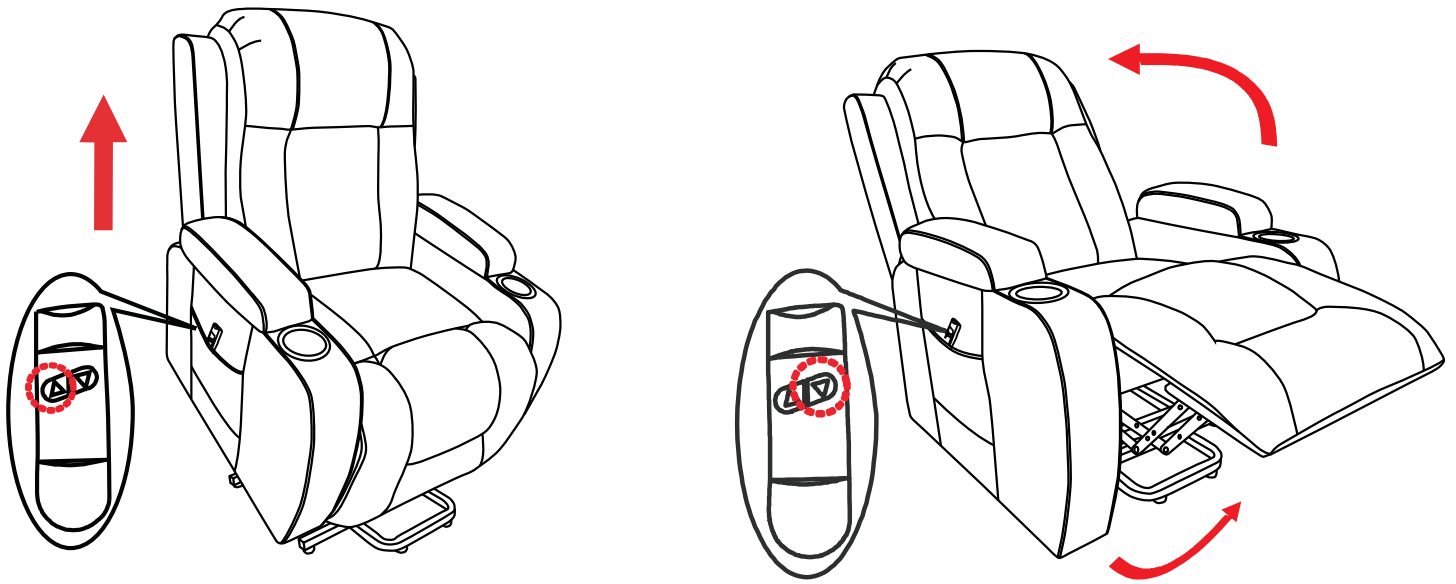Lift and Recliner Function

- Pressing the **right button**, the back will move down and the footrest will rise up. And if you press the **left button** the Lift chair will return to normal position, after that continue to press the **left button** Lift chair will rise up smoothly.
- To stop the chair at any time simply release the button tha is being pressed.

Warning

Never press both buttons together as this will damage your chair.

Never stand up while the chair is still working in order to avoiding accdent.

Troubleshooting:

If the lift chair is not working, check the following chart. Make sure you turn the power off and unplug from the wall outlet.

If you have any questions, please contact us. Do not disassemble or repair any part of the lift chair yourself.

PROBLEM

The lift chair does not function when the button on the control unit is activated.

POSSIBLE CAUSE

Check that the plug is pushed into the wall outlet properly. Check whether the power supply button is on.

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